


Professional weighing tables

For efficient and ergonomic weighing

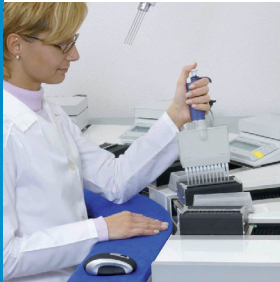

Working Ergonomically

The height of the outer table can be adjusted to individual users. The elbow cushions protect the user's arms during repetitive activities. Thanks to its insulating properties, the table is pleasant to the touch.



Low-vibration Table

The weighing table consists of two decoupled tables. The outer work surface absorbs impacts, while the balance on the inner table remains undisturbed. These versatile weighing tables are adaptable to every type of balance.



Reliable, Superior-quality Weighing

The isolation of the inner table shortens the stabilisation time of the balance. This comprehensive solution supports superior weighing performance.



First-class Design

The ergonomic weighing table saves space by providing room for a monitor and keyboard, and provides secure cable management.



Ergonomic design for every type of balance

Repetitive movements and static workplaces can create user stress. METTLER TOLEDO weighing tables are an ergonomic solution that increases weighing speed. With its adjustable height and elbow cushions, the weighing table improves work efficiency. The work surface is decoupled from the inner table providing shorter stabilisation times. METTLER TOLEDO weighing tables foster reliable weighing in your laboratory:

- Ergonomic
- Low-vibration, efficient weighing
- First-class design for professional work

METTLER TOLEDO sets the standards for ergonomic, efficient weighing.

Technical Data
Weighing Tables



	Weighing Table XP Micro	Weighing Table UMX/MX	Weighing Table XP26PC	Weighing Table MCP2
Item number	ME-11138042	ME-11138044	ME-11138041	ME-11138040
Outer table dimensions (WxDxH)	760x803x784 mm	760x803x784 mm	760x803x784 mm	760x803x784 mm
Height adjustment range (outer table)	718–820 mm	718–820 mm	718–820 mm	718–820 mm
Weight	approx. 130 kg	approx. 130 kg	approx. 130 kg	approx. 130 kg
Configuration	Outer and inner table	Outer and inner table	Outer and inner table	Outer and inner table
	Monitor base	Monitor base	Monitor base	Monitor base
	Power supply shelf	Power supply shelf	Power supply shelf	Power supply shelf
	Elbow cushion (narrow) 2x	Elbow cushion (wide)	Elbow cushion (narrow) 2x	Elbow cushion (wide)
	Stone plate	Stone plate		
	Silicone tray	Silicone tray		
			Keyboard & mouse	Keyboard & mouse

	Weighing Table MCP5
Item number	ME-11138043
Outer table dimensions (WxDxH)	1380x1380x784 mm
Height adjustment range (outer table)	718–820 mm
Weight	approx. 200 kg
Configuration	Outer and inner table
	Power supply shelf
	Elbow cushion
	Keyboard & mouse
	4-channel pump



ISO 9001 quality management certificate
 ISO 14001 environmental management certificate
 Internet: <http://www.mt.com>
 Worldwide service network

Subject to technical changes
 © 10/2007 Mettler-Toledo GmbH
 Printed in Switzerland 11795909
 MCG MarCom Greifensee

www.mt.com

For more information visit